

# WipeAlyser<sup>®</sup>

Portable Analysis Device



Available through:



## Alcohol Countermeasure Systems

60 International Boulevard  
Toronto, Ontario M9W 6J2 CANADA

+1 416 619 3500

[acs-corp.com](http://acs-corp.com)  
[drugwipe.ca](http://drugwipe.ca)


ACS, ALCOHOL COUNTERMEASURE SYSTEMS and the "Molly" are trademarks of Alcohol Countermeasure Systems (International) Inc. and are used under license.

DrugWipe and WipeAlyser are registered trademarks of Securetec Detektions-Systeme AG and are used under license.

20181002

# Ergonomic drug screening device for mobile use

It is ideal for roadside screening and drug screening in the workplace, WipeAlyser is a portable analysis device for the detection of drugs in combination with DrugWipe® S.

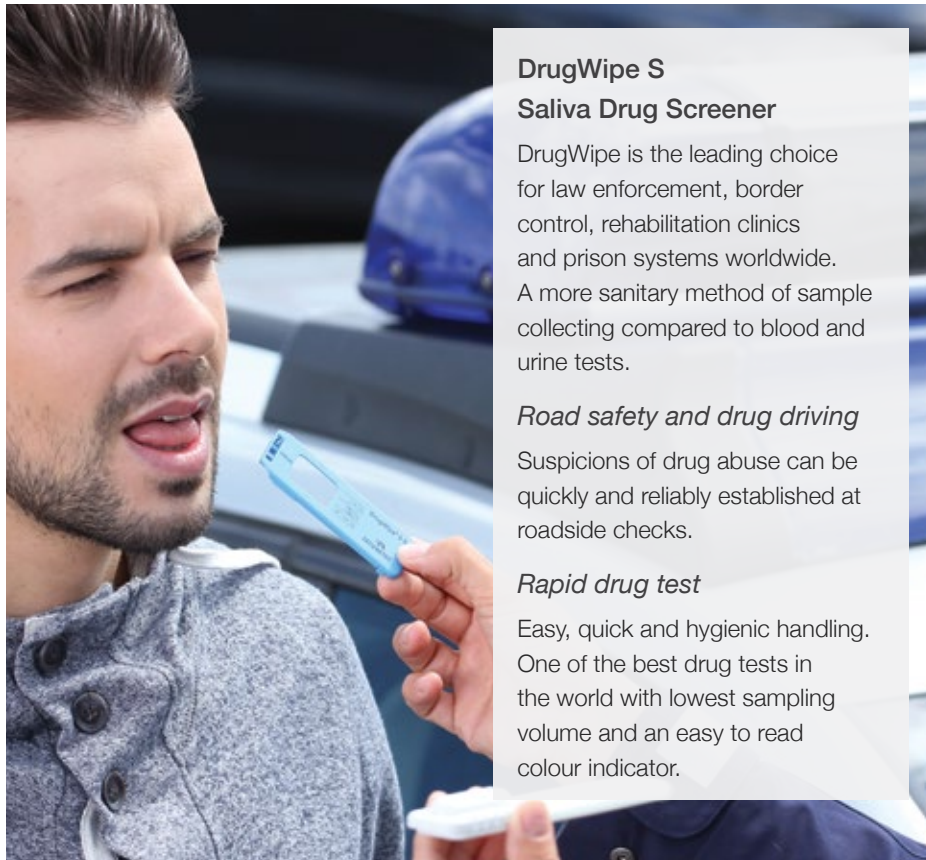

## DrugWipe S Saliva Drug Screener

DrugWipe is the leading choice for law enforcement, border control, rehabilitation clinics and prison systems worldwide. A more sanitary method of sample collecting compared to blood and urine tests.

### *Road safety and drug driving*

Suspicious of drug abuse can be quickly and reliably established at roadside checks.

### *Rapid drug test*

Easy, quick and hygienic handling. One of the best drug tests in the world with lowest sampling volume and an easy to read colour indicator.

## How does it work?

WipeAlyser is easily controlled with the integrated touch display. Entries can be made with either a finger or the supplied stylus. You will be guided through the device operation step by step. Once the prepared test cartridge is inserted into the device, analysis and reliable results are provided in seconds and read electronically.

## Why use a WipeAlyser?

### *Functional design*

- Ergonomic for mobile use
- Robust design with splash water protection
- Memory to store analysis and documentation
- Autonomy due to long battery life
- Integrated scanner for reading barcodes (optional)

### *Easy to use*

- Intuitive navigation through the analysis process
- Easy-to-operate thanks to touchscreen
- Automatic date and location via GPS



### **Display**

High-resolution colour touch display, suitable for sunlight



### **Memory**

>1,000 data records in internal memory. Additional disk space on SD card



### **Printer**

Integrated thermal printer for quick documentation



### **Ports**

USB 2.0, SD card



### **Operating temperature**

+5 °C to +40 °C



### **Battery**

Rechargeable lithium-ion



### **Battery life**

8 hours



### **Power supply**

Adapter (100-220 VAC)  
Car adapter (12 VDC)  
Rechargeable battery



### **Weight**

800 grams including printer paper roll



### **Dimensions**

24 cm x 11 cm x 8 cm