

WipeAlyser[®]

Portable Analysis Device




ACS 

ALCOHOL COUNTERMEASURE SYSTEMS

Ergonomic drug screening device for mobile use

Ideal for roadside screening and drug screening in the workplace, WipeAlyser is a portable analysis device for the detection of drugs in combination with DrugWipe® S.

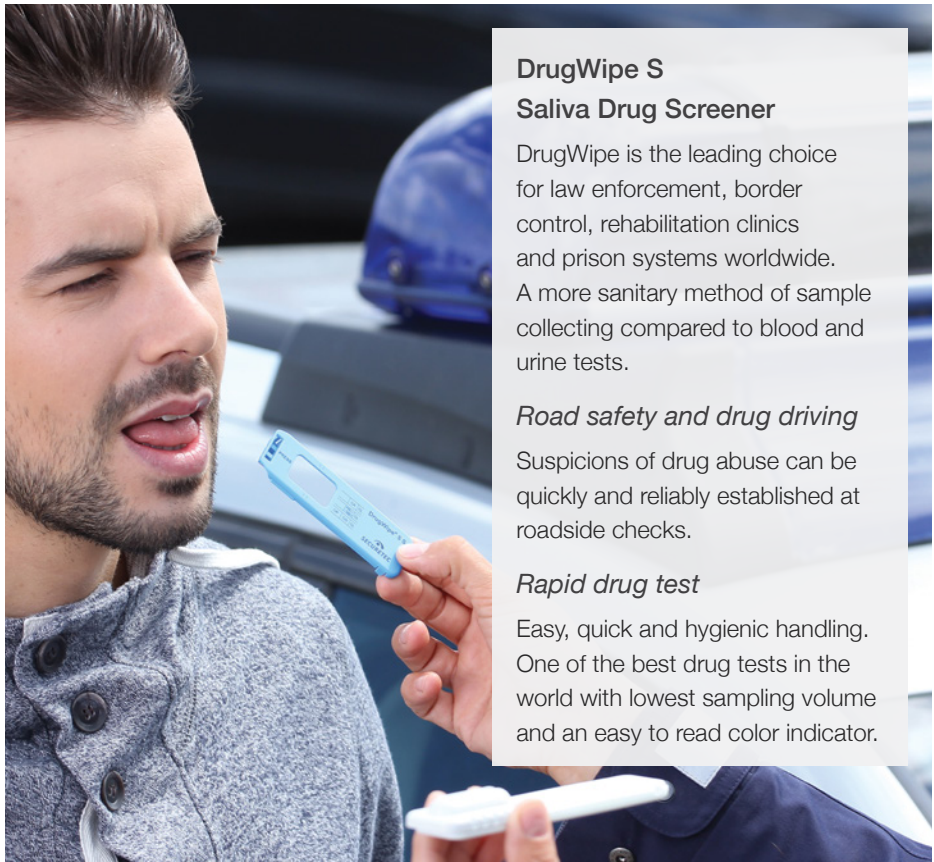

DrugWipe S Saliva Drug Screener

DrugWipe is the leading choice for law enforcement, border control, rehabilitation clinics and prison systems worldwide. A more sanitary method of sample collecting compared to blood and urine tests.

Road safety and drug driving

Suspicions of drug abuse can be quickly and reliably established at roadside checks.

Rapid drug test

Easy, quick and hygienic handling. One of the best drug tests in the world with lowest sampling volume and an easy to read color indicator.

How does it work?

WipeAlyser is easily controlled with the integrated touch display. Entries can be made with either a finger or the supplied stylus. You will be guided through the device operation step by step. Once the prepared test cartridge is inserted into the device, analysis and reliable results are provided in seconds and read electronically.

Why use a WipeAlyser?

Functional design

- Ergonomic for mobile use
- Splash-resistant and robust design
- Memory to store analysis and documentation
- Autonomy due to long battery life
- Integrated scanner for reading barcodes (optional)

Easy to use

- Intuitive navigation through the analysis process
- Easy-to-operate thanks to touchscreen
- Automatic date and location via GPS



Display

High-resolution color touch display, suitable for sunlight



Memory

>1,000 data records in internal memory. Additional disk space on SD card



Printer

Integrated thermal printer for quick documentation



Ports

USB 2.0, SD card



Operating temperature

+5 °C to +40 °C
(+41 °F to +104 °F)



Battery

Rechargeable lithium-ion



Battery life

8 hours



Power supply

Adapter (100-220 VAC)
Car adapter (12 VDC)
Rechargeable battery



Weight

28 ounces including printer paper roll



Dimensions

24 cm x 11 cm x 8 cm
(9.45 in x 4 in x 3 in)

Available through:



Alcohol Countermeasure Systems

acs-corp.com

60 International Boulevard
Toronto, Ontario M9W 6J2 CANADA
+1 416 619 3500

drugwipe.ca

5776 Hoffner Avenue, Suite 303
Orlando, Florida 32822 USA
1 877 227 7112

drugwipeusa.com

ACS, ALCOHOL COUNTERMEASURE SYSTEMS and the "Molly" are trademarks of Alcohol Countermeasure Systems (International) Inc. and are used under license.

DrugWipe and WipeAlyser are registered trademarks of Securetec Detektions-Systeme AG and are used under license.